











Quick Guide – EN

Molift Transfer Pro

molift®
by Etac

QG70101 - 2025-10-07 - Rev. 2.0

This quick guide does not replace the product's user manual. It is important to understand the contents of the user manual when operating the product. You can download the user manual at etac.com.

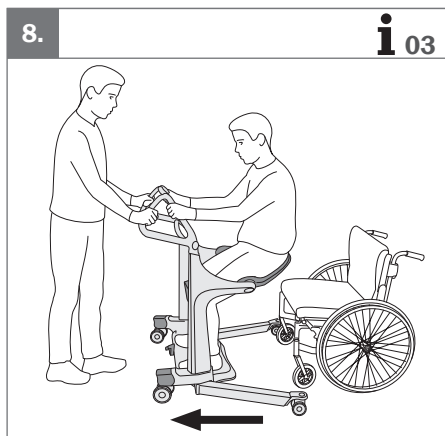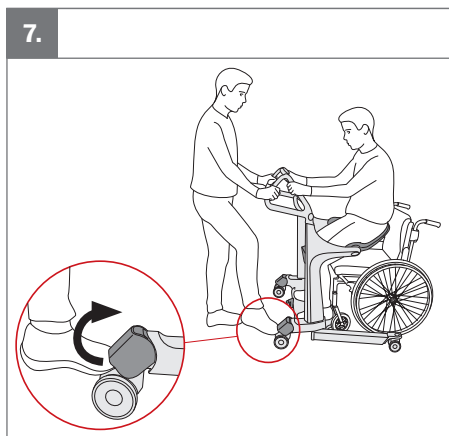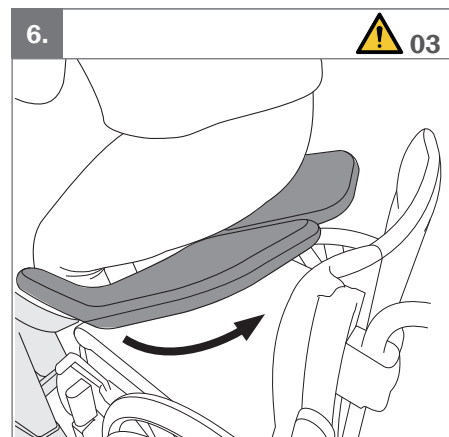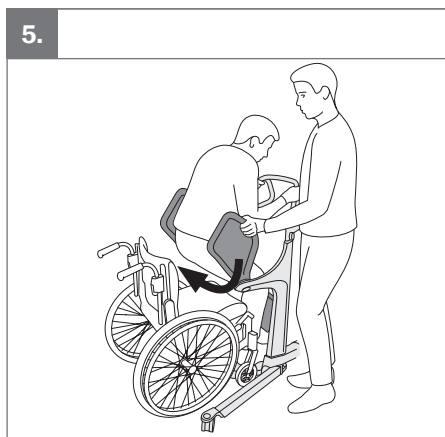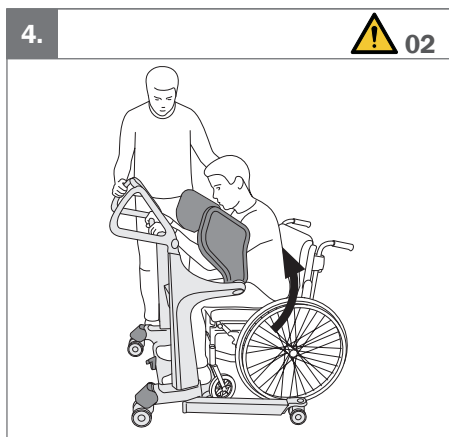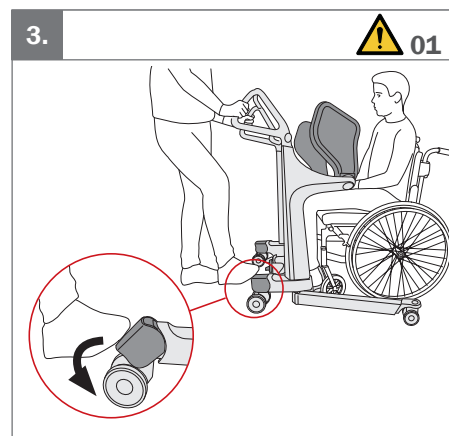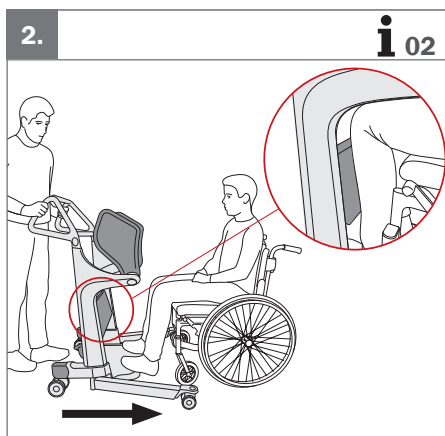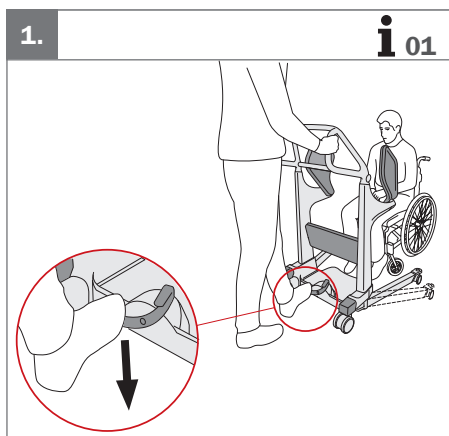










- ! Warnings**
- 01. The wheels must always be locked when the user moves between sitting and standing position.
 - 02. The user should always hold on to the handles during transfer.
 - 03. The user must remain standing until both seat pads are rotated down.

- i Information**
- 01. Adjust the legs adjustment using the pedal. Left: wider. Right: narrower.
 - 02. Place Transfer Pro so leg support reaches the user's shins.
 - 03. It is optional to move the legs in during transfer.



www.etac.com



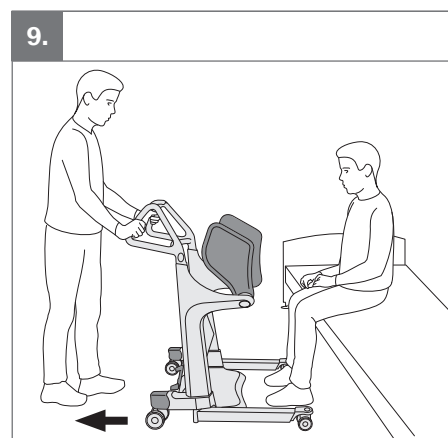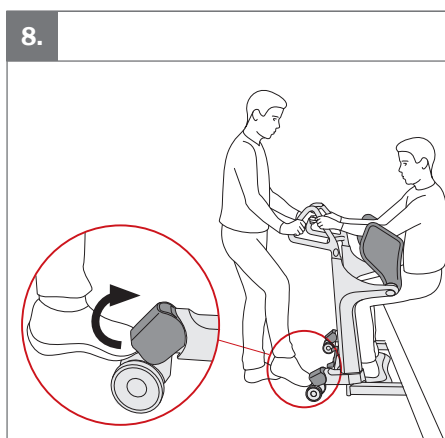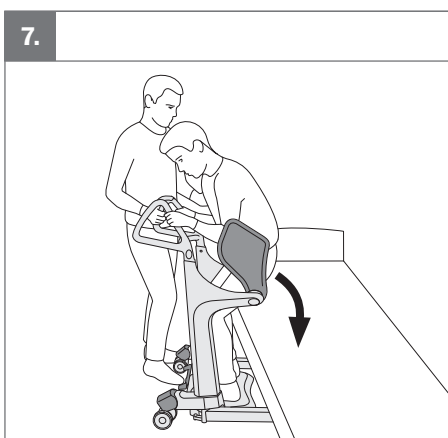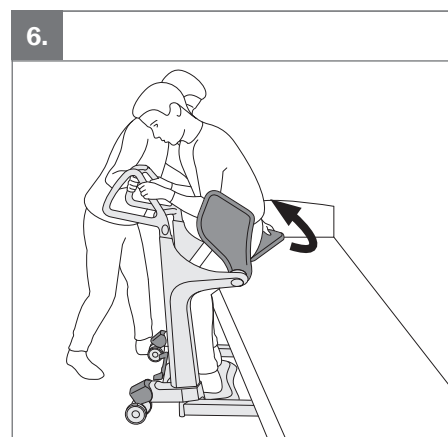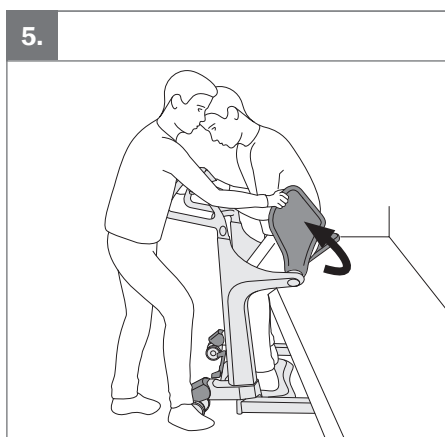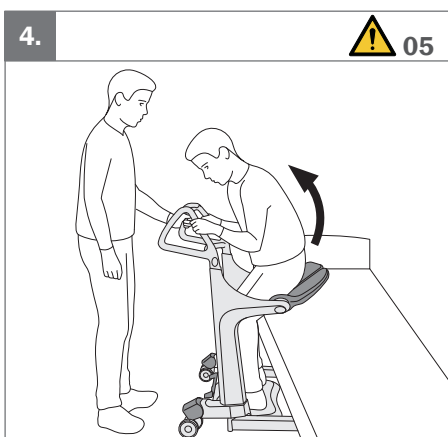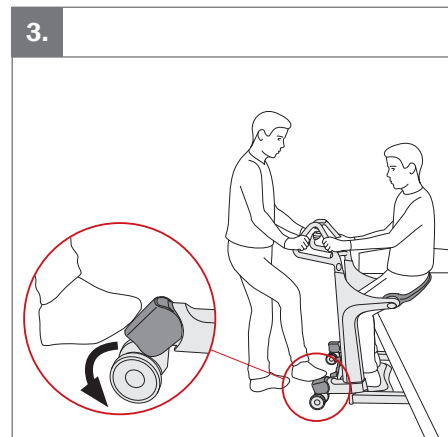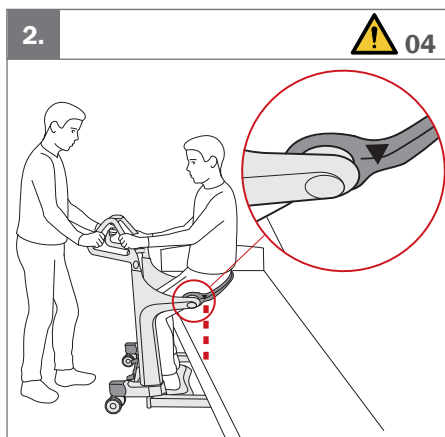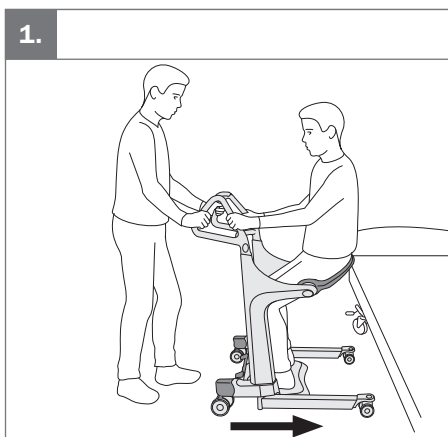

Warnings
 04. Ensure that the marking on the seat pad minimum aligns with the seat's edge before the user sits down.
 05. Always rotate both seat pads up before the user sits down.

Technical data

Weight 28,2kg / 62lbs
Max weight 200kg / 440lbs
Max height 200cm / 149cm

Cleaning

Cleaning instructions
 max 70 °C / 10 min
 pH 5-9

UK REP

Etac Ltd.
 Unit 60, Hartlebury Trading Estate,
 Hartlebury, Kidderminster,
 Worcestershire, DY10 4JB
 +44 121 561 2222

CH REP

Promefa SA
 Avenue de Morges 33
 1027 Lonay, CH
 +41 44 872 97 79



Etac A/S
 Parallevej 3
 DK-8751 Gedved

+45 796 858 33
 patienthandling@etac.com
 www.etac.com

molift®
 by Etac