



# Etac Supporter toilet arm supports

Designed to assist users and carers with every day routines.

 **etac**<sup>®</sup>

Creating Possibilities

# Etac Supporter toilet arm supports

Stable and functional toilet arm supports which are mounted directly on the toilet.

## A functional arm support

Supporter is a comfortable and stable arm support. The arm supports are foldable and have comfortable rounded grips for the hands. With its ergonomic design, Supporter blends naturally into the bathroom. The surfaces are smooth and easy to clean.

## Easy and flexible mounting

Supporter is simple to mount directly on the toilet, without the need for any tools. It fits both floor and wall mounted toilets. If space around the toilet is limited either arm support can easily be detached.

## Accessories

See our wide range of accessories on page 33 (Etac Supporter Adjustable toilet arm supports)



### Description

Supporter

### Weight

5 kg 11 lbs

### Item no.

80303112-2

### Dimensions

External width: 60 cm 23½"

Width between arm supports: 49 cm 19¼"

Arm supports height above top of porcelain: 25 cm 9¾"

Arm supports length: 41.5 cm 16¼"

### Material

Frame and arm supports: powder coated tubular steel

Arm supports cover: polypropylene

Lid and seat: melamine

Mountings: aluminium



max 85°C  
max 185°F  
for maximum of 3 minutes

max 150 kg  
max 330 lbs




Clean the product with a solvent-free cleaning agent (pH level 6–9), or with a 70% disinfectant solution. Can be decontaminated at max. 85°C / 185°F if so required.

# Etac Supporter Adjustable toilet arm supports

Adjustable arm supports that suits both adults and children.

## Flexible height and width

The arm supports can be set in three different heights and two different widths. This make them suitable for both children and adults. The arm supports can be detached, swung sideways and adjusted in height.

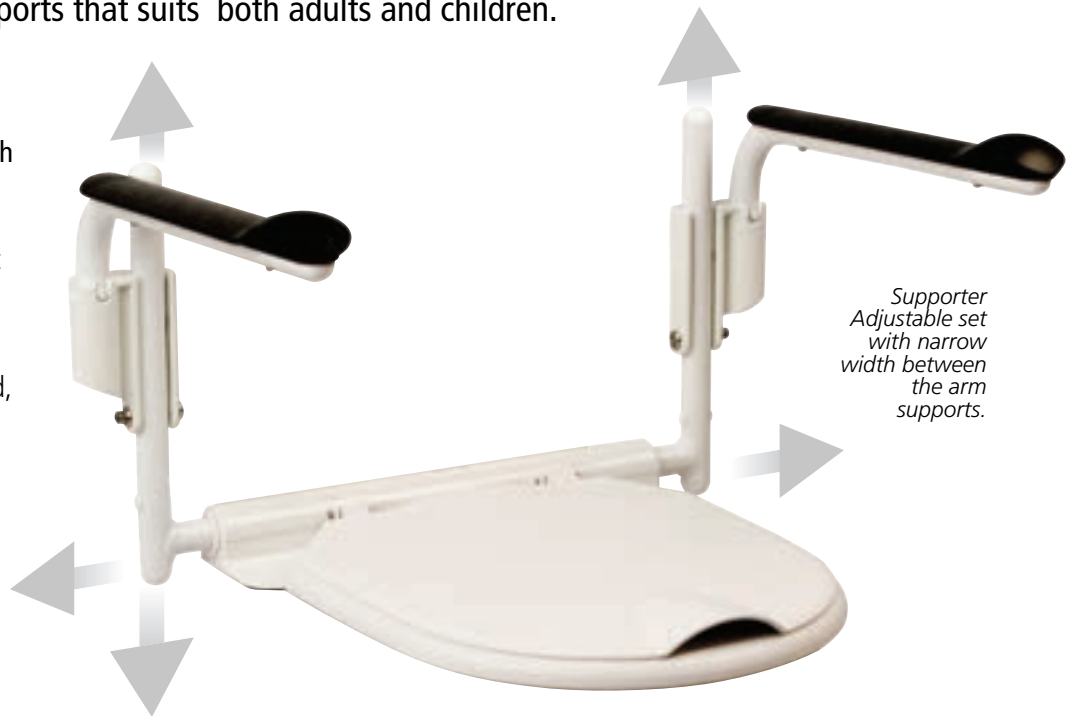

*The certified maximum weight capacity per armrest is 75kg*

## Accessories



Description	Item no.
Toilet paper holder	83030001
Compl. fitted to the arm supports.	



Description	Item no.
Crossbar 605 mm	80303034



Description	Item no.
Attachment for signal device	80303027



Description	Item no.
Attachment kit	80301052
Facilitates mounting on porcelain when space is limited for the attachment.	

Description	Weight	Item no.
Supporter Adjustable	5.8 kg 12.7 lbs	80303018-2

**Dimensions**

External width	62.3 cm 24½" to 51.3 cm 20¼"
Width between arm supports	52.5 cm 20¾" to 33.7 cm 13½"
Arm supports height above top of porcelain	20.5 cm 8" or 27.5 cm 10¾" or 34.2 cm 13½"
Arm supports length	41.5 cm 16¼"

**Material**

Frame, arm supports	powder coated tubular steel
Arm supports cover	polypropylene
Lid and seat	melamine
Mountings	aluminium


  

Clean the product with a solvent-free cleaning agent (pH level 6–9), or with a 70% disinfectant solution. Can be decontaminated at max. 85°C / 185°F if so required.



 <sup>®</sup>  
Creating Possibilities

“Everyone should be able to live a free and independent life and pursue their dreams regardless of any physical circumstances.”

 = Standard Stock Item