



novis  **etac**[®]

Etac Prio

Customisable, adaptable and scriptable

Etac Prio - Tilt in Space and Reclining Wheelchair

It's a matter of individual ability. Our aim is to make the most of the users ability, using shaping and supportive features. Accurate positioning improves the activity and mobility level. This is what we call Ability Based Seating™.

Ability Based Seating™

With age our physical abilities will decline. And abilities that are not used will decline quicker. It is important to take advantage of all possibilities to maintain optimum function and activity.

Prio is the only multi-functional wheelchair that offers a dynamic pelvic support. The support is deactivated while resting and reactivated in the active position.

The dynamic pelvic support is only one example of Prio's innovative functions that actually improve sitting tolerance over time.

We call this Ability Based Seating™.

With Prio you can shape up the position and improve quality of life.


Available in standard configuration or as a fully customised and scripted chair!



Standard Configuration Inclusions



The pelvic support is dynamic
The pelvic strap is height adjustable for precise setting. Reclining the back will automatically reduce the strap tension and retighten when the back is brought to an upright position. The ideal combination for activity and resting.

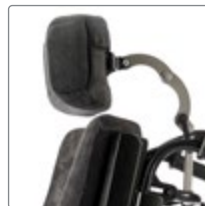

Traditional head support
The head support with turtle stem is adjustable in both height, depth and sideways. The head can be positioned in front of the trunk. This puts demands on the depth adjustment of the head support. Can be equipped with a memory function.



Attendant brake
With anti-slide handle cover. 10 cm (4") between the handle grip and brake lever.



Leg support, elevating
Swing-away and detachable. With foldable calf supports, adjustable in height and depth. Cover: Hygiene fabric.



Seat cushion, Tidy Hygiene
With hygiene cover. Material: Dartex. Water repellent, wipeable with wet cloth, washing 90° C.



Push handles
Height-adjustable and detachable, handles with anti-slide cover.

Shape Up!

All you need to know is the user's hip width. All other features can be adjusted, added, changed or replaced! The shape is created with the user in the chair. The tension adjustable back straps and back wedges are two simple techniques used to create the correct shape. Active use of the seat tilt and recline function are important.

Specific Requirements? Request a Prescription Form for customised solutions.



* Not in combination with Prio electrical tilt.

Model/Seat width	Seat depth	Total width	Transport width	Weight
40 cm (15¾")	39-51 cm (15½ - 20")	59 cm (23¼")	56 cm (22")	33.5 kg (73.9 lb)
*45 cm (17¾")	39-51 cm (15½ - 20")	64 cm (25¼")	61 cm (24")	33.8 kg (74.5 lb)
*50 cm (19¾")	39-51 cm (15½ - 20")	69 cm (27¼")	66 cm (26")	34.1 kg (75.1 lb)
55 cm (21¾")	39-51 cm (15½ - 20")	74 cm (29¼")	71 cm (28")	34.4 kg (75.8 lb)

Weights and dimensions given are for a chair with 24" rear wheels with quick-release hubs, massive tyres, stainless steel handrims, Cross brace back, leg supports, foot supports, heel straps, brakes, arm support attachment, 6½" castors with medium forks. The measurements given may vary by ± 2%.

Back support height 55 - 70 cm (21¾" - 27½") standard 45 - 60 cm (17 3/4" - 23 2/8") 3D	Balance positions Etac Prio with 20", 22" and 24" rear wheels has three alternative balance positions for each seat height.	Transport weight 22.9-23.8 kg (50.5-53 lb) without arm supports, elevating leg supports, calf supports and rear wheels.
Tilt angle -3° - +20°	Total weight Measured with arm supports, elevating leg supports, calf supports, rear wheels and anti-tippers.	Max user weight 135 kg (21.5 lb)
Back support adjustable angulation 0 - +30°	Seat height 20", 22" and 24" rear wheels: 38 - 52.5 cm (15 - 20¾") 16" rear wheels: 40 - 45 cm (15¾ - 17¾")	

*Available as a standard configuration

Options



Head support, Traditional

Also height, sideways and angle adjustable, detachable. Recommended for users with kyphosis. Cover: Grey plush or Grey hygiene fabric.



Locking pin quick release back

To enable folding the back support. Not in combination with Etac Prio electrical tilt.



Head support with side support

Adjustable in height, depth, angle and sideways. The upholstery is detachable and washable.



Seat cushion, Basic

The Basic cushion is slightly lower at the rear. Cover: Plush or grey hygiene fabric.



Memory function

Head support can be equipped with a memory function to ensure correct adjustment.



Seat cushion, Tidy Hygiene

With hygiene cover. Material: Dartex. Water repellent, wipeable with wet cloth, washing 90° C.



Back cushion

Soft and flexible back cushion, attached with velcro to be combined with side cushions and wedges. Grey plush or Grey hygiene fabric.



Seat cushion, Comfort

The comfort cushion is slightly lower at the rear. To further improve the stability, it has a removable viscoelastic foam insert and contouring for the thighs. Cover: Plush or grey hygiene fabric.



Prio 3D back support system

Unique back system providing support and space for a kyphotic back or asymmetries, attached with velcro. Can be combined with side cushions and wedges. Soft plush fabric.



Seat Velcro attachment kit

To be fastened on the solid seat and on the underside of the seat cushion.




Back support Cover

A thin cover that provides a wider back. Can be combined with side cushions and wedges. Soft plush fabric.




Push handles

Height-adjustable and detachable, handles with anti-slide cover.

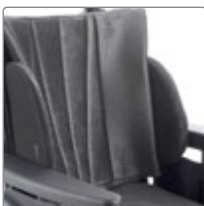


Side Cushion

Embracing support. Rotate to choose between two different shapes. Supported by the arm support. Sold per piece.

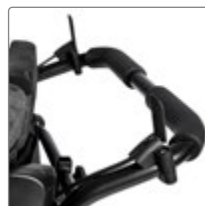

Bow handle

Height adjustable, detachable. Can be combined with attendant brake.

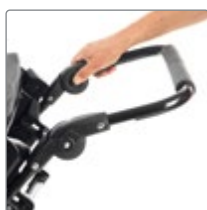

Bow handle, adjustable angle
 Height and angle adjustable. Detachable
 Can be combined with attendant brake. Fits
 Prio 40 cm, 45 cm and 50 cm seat width.



Lever locking device
 Prevents activation of the recline or tilt. To
 adjust the position, pull the locking device
 upwards, and then push it back down to lock.



Arm support
 Height and depth-adjustable,
 detachable.



Brake with angled handle
 For 20", 22" and 24" rear wheels. Provides
 a larger gripping surface. Available in black
 or red.




Arm support, narrow
 Height and depth-adjustable,
 detachable. Length: 37.5 cm (14³/₄"),
 Height: 23-32 cm (9-12¹/₂").



Brake with extension handle
 The brake lever can be angled
 forward to facilitate sideways
 transfers.



Arm support soft cover
 For added comfort on standard arm support.

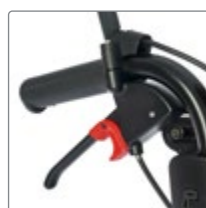

Attendant brake
 With anti-slide handle cover. 10 cm (4")
 between the handle grip and brake lever.



Arm support pad, hemi
 Adjustable in angle, the height is set with the
 arm support adjustment. Fitted onto the arm
 support narrow. Padded detachable cover in
 dark grey plush.



Upholstery sideguard
 Reduces the seat width with
 2,5 cm (1"). To be attached
 on the side guards.



Gas piston clips
 Reduce the angle options. The more clips
 added the less the chair can recline.
 No tools are required. Not in combination
 with Etac Prio electrical tilt.



**Arm support attachment
 with widening kit**
 Increases the width between
 the arm supports with 2 cm (3/4").



Comfort Lock
 This function ensures the same preset tilt
 position. There are four alternative settings
 (0°, 4°, 8° or 12°). Choose one position that
 is frequently used, for activity or rest. Not in
 combination with Etac Prio electrical tilt.



Wedges
 Cell foam with Velcro attachment. Sold in pair.
 Available in sizes Small, Medium and Large.
 Length x Width x Height
 S: 13 cm (5") x 9,5 cm (4") x 2 cm (3/4")
 M: 19 cm (7¹/₂") x 12 cm (5") x 2 cm (3/4")
 L: 20 cm (8") x 17 cm (6¹/₂") x 2.5 cm (1")

Options



Pressure distributor

Distributes pressure from the back profile.
Can also be used as a widening of the back.



Foot support, standard

Foldable, height and angle adjustable.
Equipped with heel strap, adjustable in depth.
Length: 15 cm (6").



Back angle control

Limits the back angle forward on the Prio for an active seat position.



Foot support, Cross

Height and angle adjustable. Equipped with heel strap, adjustable in depth.
Length: 13 cm (5").



Lateral support, swing-away

Adjustable in height, depth and width.
The arm support can be removed when the lateral is swung out. Ideal for side transfer.
Cushion: Polyurethane.
Cover is available as an accessory.



Foot support cover

For added comfort, plush.



Lateral support, adjustable

Adjustable in height, depth and width.
Easily removed and installed without tools.
Cushion: Polyurethane.



Calf strap

Length-adjustable and detachable.



Leg support, fixed angle

Swing-away and detachable.



Calf strap padding

Detachable. Available in lengths of 20, 30 and 40 cm.



Leg support, elevating

Swing-away and detachable. With foldable calf supports, adjustable in height and depth.
Cover: Hygiene fabric.



Padded Calf strap

Soft and warm comfort, Velcro adjustment.
Assembly on fixed or elevating leg support.



Leg support, Cross

Swing-away and detachable.



Extended foot support

Fit the Cross foot support.
Length: 21 cm (8 3/4").



Platform, foot support
The platform provides a floor sensation and stimulates movement. Foldable and detachable.



Rear wheels
20", 22" or 24". Solid tyres. Quick release hub, aluminum handrim. 16" transport wheels.



Foot support, platform Cross
Side folding, elevating. Assembly on fixed, elevating and angle adjustable leg support.



Leg abduction
Detachable and depth-adjustable.



Soft padding, platform Cross
For added comfort. Plush.



Additional seat heights for 16" rear wheels
Distance that enables two additional seat heights for 16" rear wheels: 40 cm (15¾") and 42.5 cm (16¾").




IV-stand
Comp. with attachment.



Positioning belt
Alternative fixing points on the frame.



Basket for oxygene tank
Can also be used as a carrier basket. The basket can not be combined with Etac Prio electrical tilt.

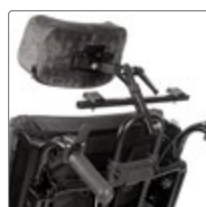

Chest strap attachment
To be combined with traditional head support.



Tray
Transparent, mounted on the arm support.



Spoke guards
Available in grey or yellow print.



Castor wheels
6" or 8".



Amputee support
Adjustable in height, depth, angle and also sideways. Fits both left and right side. Available in plush and hygiene fabric.



“Everyone should be able to live a free and independent life and pursue their dreams regardless of any physical circumstances”

novis **etac**[®]



1300 738 885 / novis@novis.com.au
12/12 Mars Rd, Lane Cove West, NSW 2066

Scan QR code for more information
novis.com.au/prio

