



Transfer Pro

Enhanced seated transfers – empowering mobility

molift[®]

Molift Transfer Pro

A transfer is more than just a movement; it is an essential part of the patient's rehabilitation and the maintenance of physical abilities. Molift Transfer Pro supports safe, dignified transfers while contributing to a healthier working environment for caregivers.

Optimised for seated transfers

Molift Transfer Pro is a sit-to-stand device designed for seated transfers between bed, wheelchair, toilet, and similar surfaces. The user should be able to follow instructions, bear some weight, and maintain core stability.

The device is suitable for use in a wide range of healthcare environments, including hospitals, nursing homes, and home care settings.

Comfort, safety, and usability

Molift Transfer Pro features an ergonomic handle with multiple grip options for both user and carer, supporting a natural movement and allowing easy access to the basin. Stable seat pads provide a secure and comfortable transfer, while the softly padded leg support adapts to the shins for increased comfort.

Easy foot placement is ensured by a low step-on height, sloped edge, and a high-friction footplate. Adjustable leg width, close positioning to the seating surface, and a one-foot-operated central brake with visual and audible feedback contribute to safe and efficient handling.

Attachment points for Molift Assist accessories are located on the side of the handle.



Handle

Multiple ergonomic grip options for the patient and carer, ensuring both comfort and ergonomic positioning of hands.



Central brake

Intuitive one-foot touch central brake, locking both wheels, regardless of their positioning, while also providing both auditory and visual feedback.



Seat pads

The stable seat pads enhance a comfortable transfer for the user. The incut edges provide a good grip, making folding them up and down easy for the carer.



Turning radius

High quality wheels for exceptionally easy manoeuvring meaning minimal force is needed when pushing/pulling the device.

Clinical benefits:

- Reduces work-related injuries or strain on the carer.
- Enables early rehabilitation.
- Comfortable and safe seated transfers.
- Supports the user to use their physical ability.
- Enables single-handed care.



Natural movement pattern

A natural movement pattern refers to the way humans naturally move during everyday activities such as standing up and sitting down. Molift Transfer Pro supports this by allowing the user to place their feet underneath the body and lean forward during the sit-to-stand movement, while the spacious handle supports a stable and ergonomic grip.



Scan the QR-code to watch a video and learn more about Molift Transfer Pro.

Product

Molift Transfer Pro

Item no.

16090300

User height

Minimum - Maximum

1490-2000 mm

Weight

Total weight

28.2 kg

Maximum user weight (SWL)

200 kg

Material

Powder coated aluminum, plastic, steel, polyuretan foam

Specifications

Height	1179 mm
Length	920 mm
Width	634 mm
Foot plate height	56 mm
Turning radius with open/closed legs	1030/1146 mm



Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution.



Molift Soft Handles 2 pcs

Helps the user reach the handle.

Molift Soft handles.....82511



Molift Assist Strap

A safety strap for extra support and safety.
Available in two sizes, S/M and L/XL.

Molift Assist Strap S/M..... 82508

Molift Assist Strap L/XL 82509




Molift Assist Sleeve

An Assist Strap with sliding sleeve, to assist user to rise up/sit down and support them during transfer.

Molift Assist Sleeve S/M..... 82528

Molift Assist Sleeve L/XL 82529

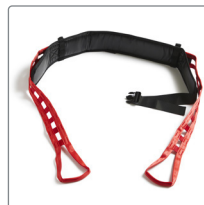


Molift Assist Belt

An assist Strap with a sliding sleeve with buckle to assist the user to rise up/sit down and support them during transfer.

Molift Assist Belt S/M..... 82964

Molift Assist Belt L/XL..... 82965

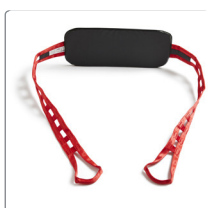

Molift Assist Sleeve Wide

An Assist Strap L/XL with firm sliding sleeve to give more support to assist user to rise up/ sit down and support them during transfer.

Molift Assist Sleeve Wide 90301



Leg support extension

Improves comfort and stability, and provides extra support for shorter users.

Leg support extension1630050



O-Bracket

The oxygen tank bracket is for securely storing and transporting oxygen tanks.

O-Bracket1630040



IV-Pole

Height-adjustable IV-Pole to be mounted on the Molift Transfer Pro.

IV-Pole1630030

Your Molift Representative

Novis Healthcare is proudly Australian Owned
For more information visit **novis.com.au**

T 1300 738 885 **E** novis@novis.com.au

20/12 Mars Road, Lane Cove West NSW 2066

molift[®]

novis[®]
enabling independence